

DSC8000

store, dry & protect your endoscope
for patient safety



Soluscope

DSC8000



a high level

Safety



*The unique system of protection
and safety for your endoscope*

AN INDIVIDUAL SHOCK PROTECTIVE TRAY

Each endoscope is securely stored in a specifically designed tray. Its profile guarantees air circulation and efficient drying.
The tray concept avoids direct endoscope handling, therefore reducing the risk of cross and environmental contamination.


YOUR SCOPES MICROBIOLOGICALLY SAFE

HEPA filtered air is flushed through all channels while air circulating around each endoscope benefits from germicidal UV-C lamp effects in each shelf.

CE MARK / EN 16442* compliant

an increased

Productivity

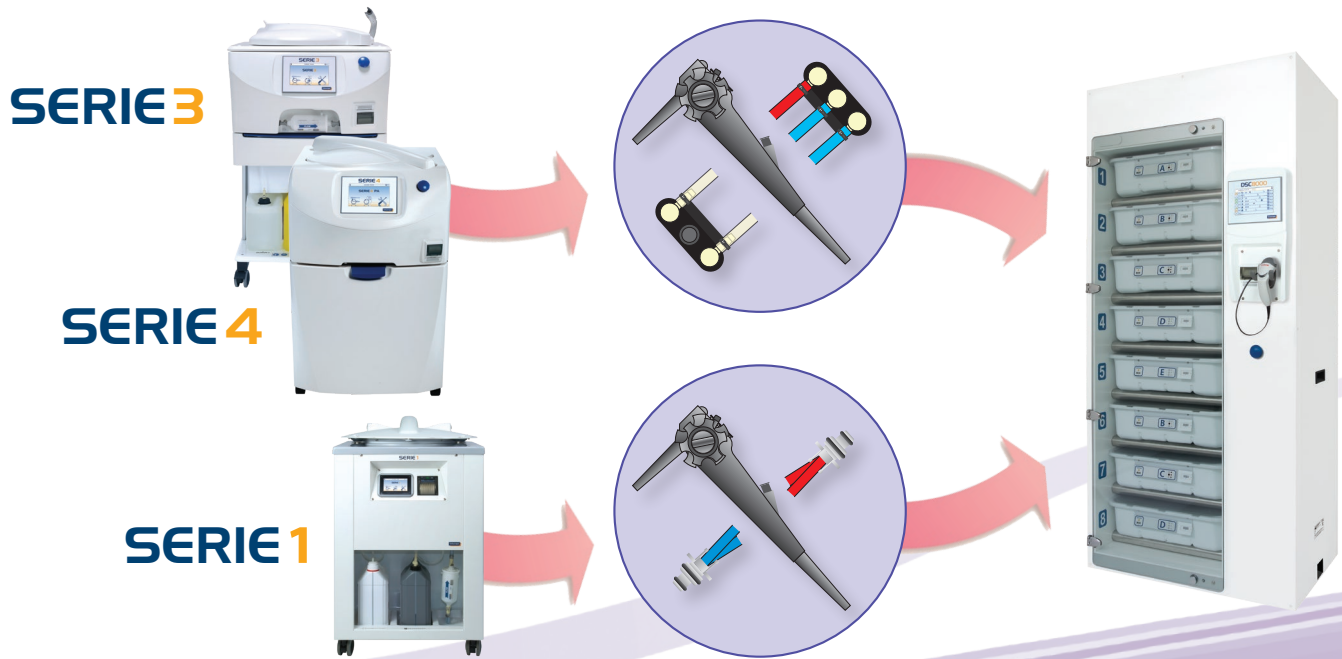

Optimize your service workload

WORKFLOW OPTIMIZATION

Your endoscope is ready-to-use in the morning, on weekends, including emergencies. The *DSC8000* interface provides real time storage data for each endoscope. You can instantly prioritize the use of your scopes and optimize your workflow.

IN-DEPTH MONITORING

The *DSC8000* provides a high level of safety by monitoring all critical parameters in real time: connections to the endoscope, air pressure in the channels, filter clogging, UV-C lamp failures.

SHARE SOLUSCOPE AER CONNECTORS TO SAVE TIME

The connector stays with the endoscope throughout the complete cycle from reprocessing to storage avoiding misconnection and reducing handling.

From *Soluscope AER* to *DSC8000*: only 1 connector.

drying and storage cabinet for endoscopes

an

Easy-to-use

device



*Offer an enhanced
comfort to your
team*

AN EFFORTLESS DEVICE

The DSC8000 brings the benefit of a new intuitive and graphic interface for fast adoption and quick use. The removable trays enable quick handling and connection to the endoscope. The large tempered glass door provides easy access and immediate identification of the scopes.

COMPLETE TRACEABILITY

The embedded traceability system includes a barcode reader for identification of scopes and operators, printer and flash memory containing electronic tickets. The DSC8000 is compatible with the IT Soluscope data management system.

A complete solution to optimize your organization

TRANSPORT

Soluscope Trolley is designed to dispatch endoscopes from reprocessing to examination room.

The specific trays equipped with transparent lids safely store the endoscopes during transportation.

Soluscope transport solution facilitates your daily routine in a safe manner.

STORAGE

Soluscope Rack can store up to 10 trays in a reduced space.

SERVICES

After-sales service, maintenance and training are designed to get the best of your device.

TECHNICAL SPECIFICATIONS

Capacity	8 endoscopes
Validated storage time	Up to 7 days
Endoscope registration	Up to 150
Electrical supply	220-240V / 50 Hz or 220V / 60 Hz
Air supply	Embedded compressor or connection to medical / technical air
Air treatment exterior endoscope surface	HEPA Filter & Germicidal UV-C lamps
Air treatment in the channels	HEPA Filter (compressor version) or Medical / technical air
Traceability	Barcode reader - Printer - Removable memory card - IT Soluscope (option)
Version	Single door / Pass-Through
Dimensions (W x D x H)	90 x 64 x 202 cm (compressor) 90 x 64 x 205 cm (medical / technical air)
Weight	170 kg

Data management solution

The *DSC8000* is compatible with the global traceability software solution *IT Soluscope* which provides access to your endoscope data base in real time allowing data management.



Soluscope
REGISTERED TO ISO
9001:2015 AND 13485:2016
FILE NUMBER A14025



SOLUSCOPE COMPLETE SOLUTION

Soluscope proposes a complete solution for flexible endoscope reprocessing.

It offers products for every step of the reprocessing cycle.

Each part is certified, easy to use and effective in order to reduce your daily load.



www.soluscope.com

Soluscope SAS : 100, rue du Fauge - Z.I. Les Paluds - 13400 Aubagne - France

☎ +33 (0)4 91 83 21 22 📠 +33 (0)4 91 83 21 10 ✉ info@soluscope.com

Soluscope products are compliant with *European and international standards*, and are *CE marked*. *ISO 9001 version 2015* and *ISO 13485 version 2016* certified, *Soluscope* is continuously committed to satisfying its clients rigorous demands.

distributed by